

Dual Convection Oven YEAR 2023

DOF-H45 II DOF-H80 II DOF-H140 II DOF-H250 II



Microprocessor controller with LCD display, accurate and reliable.



High quality stainless steel chamber, removable shelf, easy-to-clean.



Silicon sealing ring for reliable sealing.



Equipped with leakage protection.




Anti-hot handle



Equipped with spare temperature control which ensures the product work normally even the main temp.control failed.

Features

- Equipped with multiple sets of heaters.
- Adopt temperature selection switch which can select the grade according to the heating speed and operating temperature.
- Convenient operation for air inlet, reasonable air duct structure, good temperature uniformity.

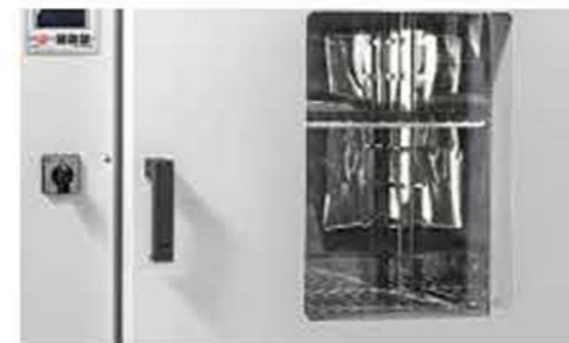

With inner glass door for easy observation.

Options

- Multi-segment programmable control
- Built-in Printer
- RS485 interface

Specifications

Model	DOF-H45 II	DOF-H80 II	DOF-H140 II	DOF-H250 II
Cycle Mode	Dual Convection: Force Air and Natural convection			
Capacity	43L	81L	138L	252L
Temp. Range	RT+10-300℃			
Temp. Fluctuation	±1℃			
Temp. Resolution	0.1℃			
Temp. Uniformity	≤Max. Temp. ±3.5℃			
Shelves	2pcs/adjustable			
Consumption	1300W	2000W	2500W	3000W
Power Supply	AC110/220V±10%, 50/60Hz			
Internal Size(W*D*H)	350*350*350mm	450*400*450mm	500*500*550mm	600*600*700mm
External Size(W*D*H)	640*500*600mm	730*530*690mm	780*690*800mm	880*730*950mm
Package Size(W*D*H)	750*630*770mm	850*660*860mm	900*760*980mm	1000*860*1100mm
Net/Gross Weight	35/55kg	50/80kg	75/110kg	110/150kg